


Cemento-MP

TUNED SYSTEM FOR VERTEBROPLASTY

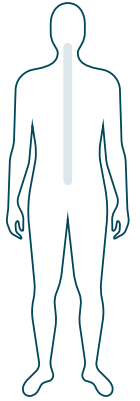

Cemento-MP

THE REUSABLE VERTEBROPLASTY SYSTEM

The vertebroplasty system consists of the **reusable application system Cemento-MP**, the special vertebroplasty needles, injection cylinder, and connection system.

Connection tube

FEATURES


Large screw applicator for a powerful and controlled injection


Special vertebroplasty needle featuring shock-resistant bayonet Luer lock fittings

Pressure-resistant injection cylinder with special reinforcement and rotating Luer lock

Reusable application system

The special vertebroplasty needle has indentations on the needle and stylet that indicate the opening of the bevel edge at the needle tip. The special bevel edge makes it easier to control the needle and the cement delivery.

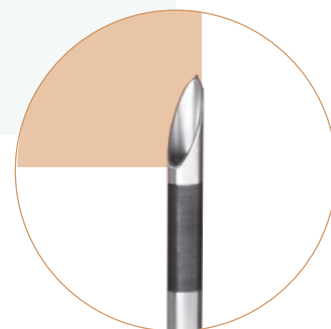

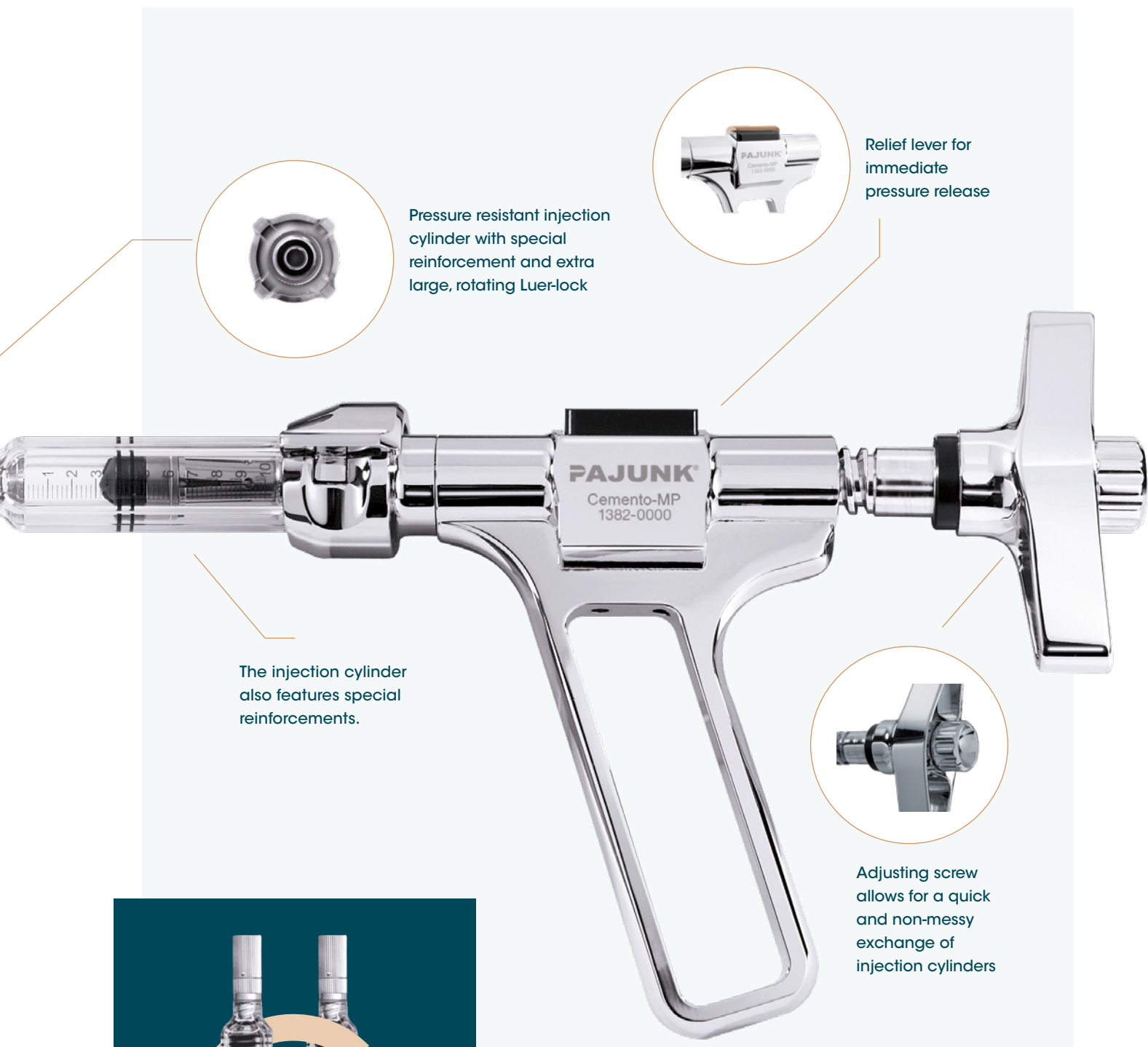

Pressure resistant injection cylinder with special reinforcement and extra large, rotating Luer-lock



Relief lever for immediate pressure release



The injection cylinder also features special reinforcements.



Adjusting screw allows for a quick and non-messy exchange of injection cylinders



RAPID EXCHANGE CYLINDER MECHANISM




Reduction adapter - Standard Luer-Lock





Aspiration cannula

CEMENTO-MP REUSABLE VERTEBROPLASTY SYSTEM

 Cemento-MP	
Order code	
Reusable application system incl. flushing adapter	1382-0000

 Cemento Mini Set for Cemento-MP	
Order code	
Supplementary set suitable for the Cemento-MP system	1382-0100

 Special vertebroplasty needle set for Cemento-MP		
Diameter	Needle length	Order code
mm / gauge	mm	
3,2 / 10	100	1382-1010
	150	1382-1510
2,4 / 13	100	1382-1013
	150	1382-1513

 Special vertebroplasty needle for Cemento-MP		
Diameter	Needle length	Order code
mm / gauge	mm	
3,2 / 10	100	1394-1010
	150	1394-1510
2,4 / 13	100	1394-1013
	150	1394-1513
1,8 / 15	100	1394-1015
	150	1394-1515

FOR THE VISCOUS BONE CEMENT CEMENTO FIXX-M* VISIT OUR WEBSITE

*Manufacturer: G-21 S.r.l., Via S. Pertini 8, 41039 San Possidonio (MO), Italy 



Order Code 9997-0019EN
Rev (07/2024) | @optimed



Distributed by:
optimed Medizinische Instrumente GmbH
Ferdinand-Porsche-Straße 11 | 76275 Ettlingen Germany
Tel. +49 (0) 7243 / 76 33-0 | Fax +49 (0) 7243 / 76 33-624
info@optimed.com | www.optimed.com

Manufactured by:
PAJUNK® Medical Produkte GmbH
Karl-Hall-Straße 1
78187 Geisingen | Germany

