
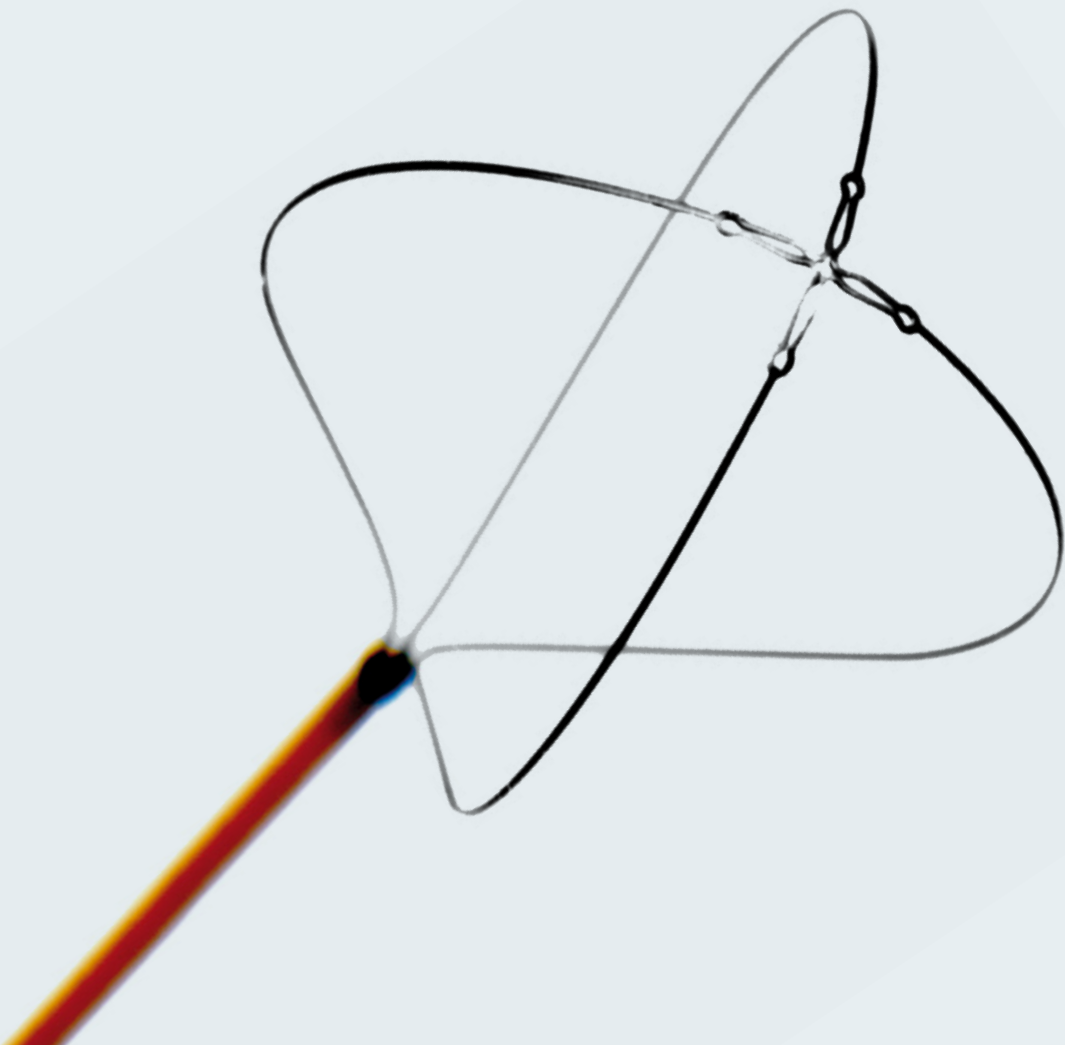


# Stone Baskets

FREEMAN, EQUADUS & DOURO

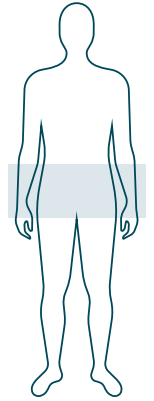


# STONE BASKETS FOR ENDOUROLOGY

FOCUSED ON THE REQUIREMENTS  
OF STONE MANAGEMENT PROCEDURES

The stone baskets were **designed for the removal of ureteral and kidney stones**. The focus was particularly on the special anatomical and interventional requirements of endourology. The unique properties of Nitinol enable shape retention with high radial force to maintain **stability and control**. The handle is designed **for precise, single-hand operation**.

## FREEMAN

DESIGNED FOR STONE  
RETRIEVAL IN THE RENAL PELVIS

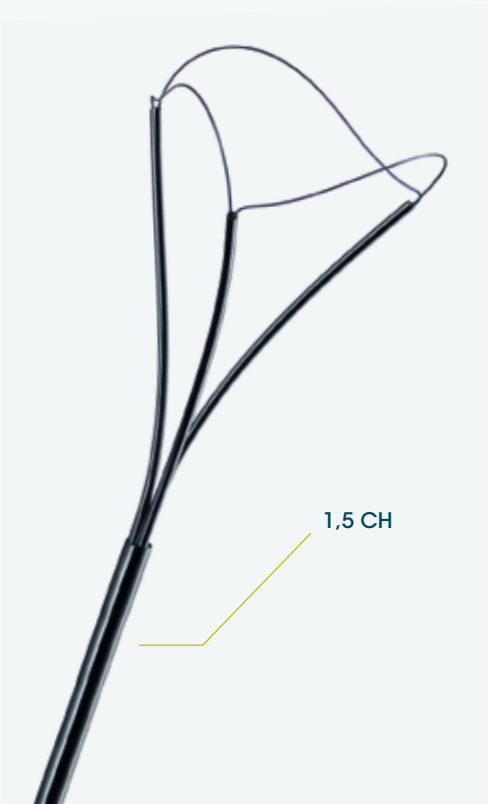


### FEATURES

Clever design: The hybrid basket offers the benefits of a grasper and a basket


Small diameter enables a high irrigation flow for brilliant visibility

Special feature: Stepless locking option of the basket on the handle

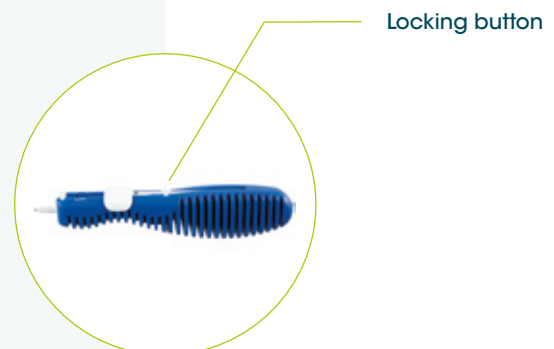

# EQUADUS

THE ALL-ROUNDER  
IN THE RENAL PELVIS

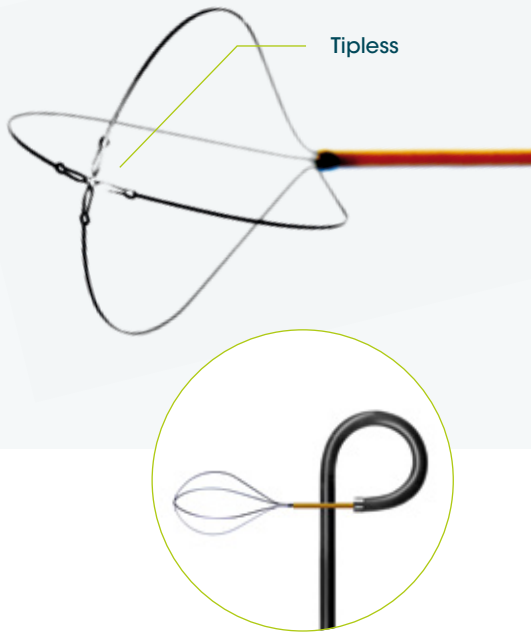

## FEATURES

Opens easily and remains torsionally stable – even when the ureterscope (URS) is flexed

Truly tipless due to the laser-cut tip

High nitinol radial force – in particular for the narrow renal pelvis

# DOURO


THE ALL-ROUNDER  
IN THE URETER

## FEATURES

High pushability – for challenging anatomies or the passage of ureteral stones

High radial force – the stone basket also opens in the narrow ureter

Basket made from nitinol – for high shape retention

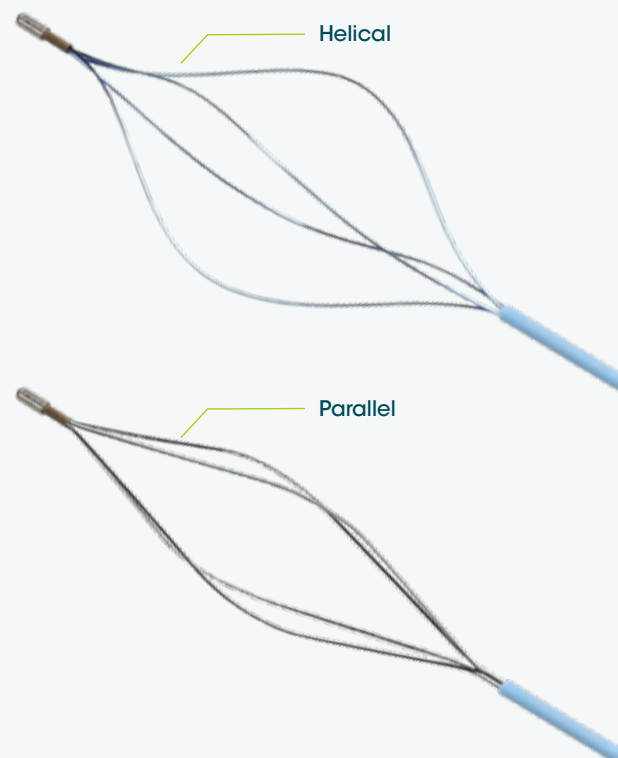

## SCOPE OF DELIVERY


- Stone basket
- Insertion aid (Freeman and Equadus)
- Tuohy-Borst adapter (Freeman and Equadus)



Packing unit:  
1 piece


| Product   | Sheath |             | Basket |             | Order Code |
|---|--------|-------------|--------|-------------|------------|
|   | Ø / CH | Length / cm | Ø / mm | Length / mm |            |
| Freeman<br>          | 1,5    | 120         | 11     | 13          | 49012690   |
| Equadus<br>        | 1,8    | 120         | 11     | 15          | 47012172   |
|   | 2,2    | 120         | 12     | 17          | 42013306   |
| Douro helical<br>  | 2,5    | 90          | 15     | 36          | 42013307   |
|   | 3,0    | 90          | 15     | 36          | 42013309   |
| Douro parallel<br> | 2,5    | 90          | 15     | 36          | 42013308   |
|   | 3,0    | 90          | 15     | 36          | 42013310   |

**FOR MORE ENDO-  
UROLOGICAL HIGHLIGHTS  
VISIT OUR WEBSITE**

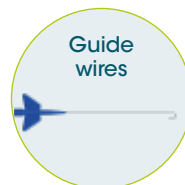

Order Code 9997-0003EN  
Rev (06/2024) | ©optimed

Distributed by:  
optimed Medizinische Instrumente GmbH  
Ferdinand-Porsche-Straße 11 | 76275 Ettlingen Germany  
Tel. +49 (0) 7243 / 76 33-0 | Fax +49 (0) 7243 / 76 33-624  
info@optimed.com | www.optimed.com

Manufactured by:  
EPflex Feinwerktechnik GmbH   
Im Schwöllbogen 24  
D-72581 Dettingen an der Erms | Germany